smartLAB[®]

walk B

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smartLAB[®]

walk B

3-D Step-counter with wireless data transfer

User Manual

Please read this user manual thoroughly before first using this device

HMM Diagnostics GmbH
D-69221 Dossenheim, Germany
www.hmm.info

Content

I. Introduction

II. Your smartLAB[®] *walk* B

Display & functions

Specifications

III. Setup & Functions

Inserting the battery

Mounting the carrying cord

Important information connecting to the fitmeFit move App and hLine App

Before the first use:

Setting weight, step length, signal tone and time:

Manually transfer Not Sent data

Manually transfer all data

Starting up

Memory function

Additional functions

Maintenance

Regularity and Safety Notice

Warranty

28

29

29

29

30

30

31

32

33

33

40

41

43

43

44

45

46

48

I. Introduction

Thank you for using the **smartLAB®walk B** 3-D step counter. Because of its small and handy form and its appealing design, this device suits best your needs for the control of walked distances. Because of the integrated 3-D sensor technology, this device counts normal steps as well as lunges or backward steps. You can comfortably put it in your pocket and it will count all your movements.


Besides, this step counter has an integrated Bluetooth Smart module, making it possible to transfer the measured data (steps, calories, distance) to an accordant receiving device like PC, laptop, tablet or smartphone.

smartLAB®walk B step counter helps you find your own individual training level and thus, enables you to control and enhance progress in your personal fitness workout.

smartLAB®walk B step counter records the number of steps per minute. This feature enables the receiving device to analyse the data received. The software can show you when and how intensive you made your steps. This can only be seen on the software of the receiving device.

II. Your smartLAB®walk B

Display & functions

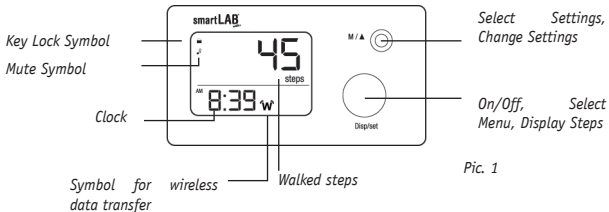
Specifications

1. Name: **smartLAB®walk B**
2. Size: 70 x 37 x 10,5 mm
3. Weight: 30 g
4. Battery: 1 x CR 2032 battery
5. Unit for steps walked: amount
6. Unit for calories: kcal
7. Unit for distance: km
8. Unit weight: kg / lb (adjustable)
9. Communication / data transfer: BLE

III. Setup & Functions

Inserting the battery

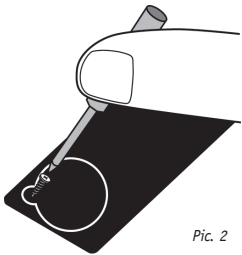Pic. 2

In order to change/ insert the battery, you need to unscrew the battery lid on the back of the device. Place the coin cell with the plus side to top into the opening. If you change the battery, use the screwdriver to lever off the battery. Afterwards, put the lid back onto the battery and screw it down again.

Note: After a battery change, you need to set the time of the device again. Your saved data will stay on your device when changing the batteries. An automatic setting of your data is enabled by certain receiving devices. Please find more information under *Heart Beat Mode*.


Mounting the carrying cord

You can wear your **smartLAB[®]walk B** with the carrying cord included in the set. To do so, you will have to connect the cord to the device by doing the following steps:

Push the thin end of the cord through the two tiny eyes at the back side of the device. Pushing together the cord a little bit makes it easier to squeeze it through the eyes (pic. 3). After that, push the thicker end of the cord through the loop on the thin end and pull it tight (pic. 4)

Pic. 3

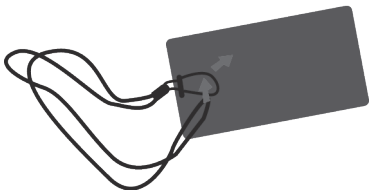Pic. 4

Important information connecting to the fitmeFit move App and hLine App

You can transfer your data to the fitmeFit move App or hLine App, available for Android and iOS. Please note following information:

- The connect of the smartLAB walk B to your mobile device is possible through the App and not through the Bluetooth settings of your device
- To transfer data the serial number of the device needs to be entered into the App, either automatically or manually
- The setting information can be automatically send to your step counter thanks to the Heart Beat Mode, e.g. weight or step length. Manually entered information will be then overwritten
- With iOS a data transmission is only possible when the App is activ and in the foreground
- You can see in the App your daily walked steps as graph with more detailed information. Hold your mobile device transversely and then choose a day which you want to see in details

Before the first use:

Before first using your **smartLAB®walk** B step counter, please hold the „Disp/Set“ button pressed for more than 3 seconds. The device will go on and enters the setting mode.

The **smartLAB®walk** B step counter has an automatic setup function (Heart Beat). This way there is no need to setup the step counter manually. Please refer to Heart Beat section for more details.

Setting weight, step length, signal tone und time:

Setting weight:

1. Press the “M/△” button and choose between „kg“ and „lb“. In order to confirm your choice, please press the „Disp/ Set“ button.
2. After you choose your unit, the weight indicator will flash. Now you can enter your personal body weight and confirm by pressing the “M/△” button. The possible weight range is between 30 and 200 kg.

With the “M/△” button you can adjust your weight. You can only increase the value. After reaching the maximum level of 200 kg, **smartLAB[®]walk B** will start with 30 kg again.

By keeping the “M/△” button pushed, the weight value increases faster. After having set your personal weight, press „Disp/ Set“ to confirm.

Setting step length:

After you entered your body weight, the step length symbol flashes. Now you can enter your step length by pressing the “M/△” button. The step length can vary between 30 and 120 cm. After having set your correct step length, confirm with „Disp/ Set“.

Information on determining the correct step length:

To determine the correct the correct step length please follow these steps:
Choose a certain distance you know the length off. Walk the distance and count your steps (e.g. with the step counter). Now calculate the step length the following way: the length of the walked distance (in cm) divided through the number of walked steps.

Setting signal tone:

After you entered your step length, you can set the device signal tone (🔊). By pushing the “M/△” button, you can choose if the signal beep is activated (“ON”) or not (“OFF”). Please confirm again with the „Disp/ Set“ button.

Setting time:

After that, your device enters the time setting mode. Press the “M/△” button in order to switch between the two display modes 12h and 24h (24h mode is preset).

By pushing „Disp/ Set“ button you can set the time. First you set the hours by pushing the “M/△” button. Confirm by pushing „Disp/ Set“. Afterwards, you can set the minutes by repeating the steps mentioned. Finally confirm again with the „Disp/ Set“ button.

Note: Setting your weight and step length is essential for the display of walked distances and burned calories. You need to enter correct values for both parameters. Otherwise it is impossible for **smartLAB[®]walk B** to display correctly your walked distance and calory consumption.

Setting date:

After that, your device enters the date setting mode. Press the “M/△” button to enter the current day (display shows “D”) and confirm with “Disp/ Set”. Proceed accordingly for setting the month (display shows “M”) and year (display shows “Y”). Always confirm with “Disp/ Set”.

Note: In order to use wireless data transfer, date and time need to be set correctly. Otherwise, measured data can not be assigned accordingly.

Bluetooth Smart settings:

Bluetooth Smart is a wireless data transfer standard. In order to use it, transmitting- and receiving device must be equipped with a specific module. Your **smartLAB[®]walk B** step counter has an integrated Bluetooth Smart module, making it possible to transfer measured data to a PC, laptop, tablet or a mobile phone. You can store, analyse and process the data with the accordant telemedicine software on these receiving devices.

You have 2 different Bluetooth Smart setting modes:

1. Per 10 minutes mode (recommended)
2. Heart Beat mode

Setting the Per 10 minutes BLE mode:

In this mode data will be sent out every 10 minutes to the accordant receiving device. All data which was stored during the last time interval will be sent out then. You can see the automatic mode in picture 5. Switch the mode ON or OFF by using the “M/△” button.

It is recommended (Factory Setting) to keep this mode (ON) to enable data transfer.

pict. 5

Setting the Heart Beat mode:

The Heart Beat mode allows an automatic setup setting such as date, time, weight and step length. These information are provided by the receiving device/ software. It is important that the receiving device supports the heart beat mode. If entering your profile information e.g. into the fitmeFit move App or hLine App, these information will then be transmitted to your **smartLAB** walk B step counter.

If this mode is (ON), which is the factory setting, the step counter will try to communicate every minute with a receiving device. Is there any receiving device in reach which supports this function, the step counter will be set automatically. After the step counter was set, the heart beat function will start again every 30 minutes.

The benefit of this function is to update the setting automatically. Another benefit is that the step counter sends the battery status to the receiving device. This way the receiving device can warn the user that the battery is almost empty. You can see Heart Beat in picture 6. Switch the mode (ON) or (OFF) by using the “M/△” button.


pict. 6

Manually transfer Not Sent data

You can also manually trigger data transfer to the receiving device. To do so, you first have to turn on SEND mode on your **smartLAB[®]walk B** device.

Please quickly press and release the “Disp/ Set” button, and then press it again for 6 seconds. In the upper display area “SEnd” appears.

If you now press „M/△” the display will show “picture 7”.

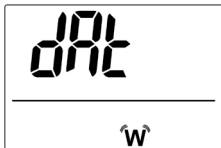pict. 7


By pressing “Disp/ Set” again, you will manually trigger the data, which has not been sent yet, transfer. The symbol for wireless data transfer **W** will flash in the display. If the data transfer was successful, you will see “picture 9”; if the data transfer could not be completed, “picture 10” will appear on display. In case no connection at all can be established, the display will show “picture 11”.

Manually transfer all data


You can also manually trigger all data transfer to the receiving device. To do so, you first have to turn on SEND mode on your **smartLAB[®]walk B** device.

Please quickly press and release the “Disp/ Set” button, and then press it again for 6 seconds. In the upper display area “SEnd” appears.


If you now press „M/△” twice the display will show “picture 8”.

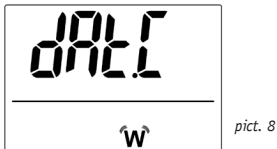
By pressing “Disp/ Set” again, you will manually trigger the all data transfer. The symbol for wireless data transfer **W** will flash in the display. If the data transfer was successful, you will see “picture 09”; if the data transfer could not be completed, “picture 10” will appear on display. In case no connection at all can be established, the display will show “picture 11”.

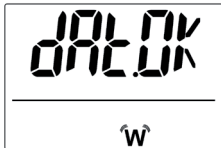pict. 09

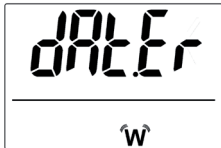pict. 10

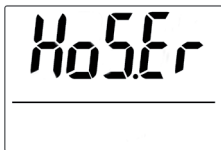pict. 11

Starting up

- After having set time and step length, you can use your **smartLAB[®]walk B** step counter. Just take the device and start walking. The device will register a movement.
- After the 10th step, the counter starts counting your steps.
- To see how many steps you have walked so far you need to turn on the Display. Press on the „Disp/Set“ button to unlock the Display. It will then turn on.
- After physical exercise or long walks, simply press „Disp/ Set“ to view your total step count, walked distance, time and total calory consumption.

Memory function

- Your **smartLAB[®]walk B** step counter automatically saves your records on a 35-day basis.
- You can see the stored data by pressing the „M/ Δ “ button for a few seconds. Please choose the desired day in the according display. You can now see a total

step count, total distance, the time of your exercise and the calory consumption.

Note: Once you have walked 99999 steps, the device stops counting your steps. The LCD display will start twinkling and reset after the time setting was set back to 0:00.

Additional functions

Battery-display:

When the low battery symbol  appears in the display, your **smartLAB[®]walk B** needs new batteries. Please exchange the batteries.

Step counter and batteries need to be disposed of according to local regulations.

Exchanging the battery:

When exchanging the battery, the display will blink for 5 seconds. After that you can use it normally.

Locking your device:

If no button is pushed for 30 seconds, your **smartLAB®**walk B automatically activates key lock. This is shown by the according display symbol. By pushing the „Disp/ Set“ button, you can unlock the device.

Display-power-off:

When you are in user mode, your device automatically turns off after 1 minutes of non-use. To turn it back on, click on the „Disp/Set“ button.

Maintenance

- Do not dump the device or use force while operating it
- Avoid exposing the device to high temperatures or direct sunlight.
- Do not immerse in water, since this might cause irreparable damage.
- Let the device acclimate after having stored it in a temperature close to zero. Do not immediately use it in such a condition.
- Do not try to disassemble the device.
- If you do not use the device for a longer time, please take out the batteries.

Regularity and Safety Notice

Declaration of Conformity

This product is compliant with the EC Directives: -1999/5/EC

Safety Statement

This product has been verified to conform to the safety certification in accordance with the specification of EN Standards: EN 60950-1:2006 + A11:2009 + A1:2010 + A12:2011

WEEE note

The WEEE (Waste Electrical and Electronic Equipment) Directive, that came into effect as European law on 13th February 2003, led to a major change in the disposal of electrical equipment.

The primary purpose of this directive is to prevent of electronic waste (WEEE), while encouraging for reuse, recycling and other forms of reprocessing to reduce waste.

The WEEE logo on the product and the packaging indicates that it is not allowed to dispose of the product in the ordinary household waste. It is your responsibility to deliver all the disused electrical and electronic equipment to the respective collection points. A separate collection and proper recycling of electronic waste helps dealing economically with natural resources. Furthermore, the recycling of electrical waste is a contribution to the conservation of our environment and thus the health of all people. More information concerning the disposal of electrical and electronic devices, reprocessing and the collection points you can get from the local authorities, waste disposal companies, from retailers and manufacturers of the device.

RoHS compliance

This product complies with Directive 2002/95/EC of the European Parliament and the Council of 27th January 2003 regarding the limited use of dangerous substances in electrical and electronic equipment (RoHS) and its variations.

Warranty

HMM Diagnostics GmbH products need to fulfill high quality requirements. Because of this reason, HMM Diagnostics GmbH gives a 2-year warranty by purchasing this **smartLAB**[®] product. You can even extend the warranty from 3 to 5 years without extra pay when you register your product. Please register online under following website: www.hmm.info/en/registration.

Wear parts, batteries etc. are excluded from warranty.

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